

SETTING UP the Eversense 365 App



1 Create your Eversense account


- Download this app → 

You can download the free **Eversense 365 App** from the Apple App Store or on Google Play*.

- When you open the app, accept the License Agreement and tap **Create an account** on the LOG IN screen.
- Complete the information on the **Create an account** screens
- When you are finished, click **SUBMIT**.

2 Activate your Eversense account

Once you have created your account, you will receive an email from Eversense with a link to complete the registration and activate your account. You can complete this process entirely on your mobile phone.

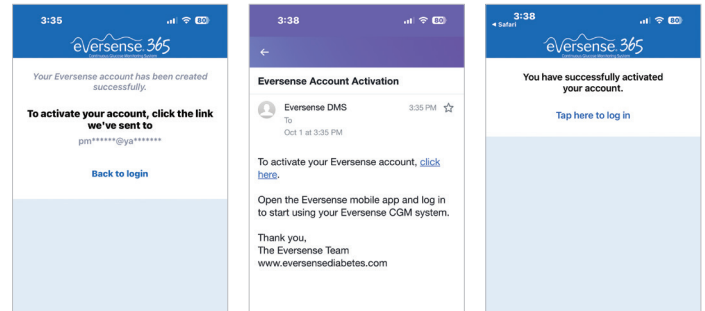

Important: You must activate your account within 72 hours by clicking on the link in the email.

Congratulations! Your account is now activated.

3 For future reference, write down your Eversense account log-in information below:

Email: _____

Password: _____

Scan this QR Code to watch a video on how to set up the Eversense 365 App



To contact Customer Support, call 844-SENSE4U (736-7348)

* For a full list of compatible devices, please visit <https://www.eversensecgm.com/compatibility/>

Senseonics

© 2026 Senseonics, Inc.
All rights reserved.

Manufactured by Senseonics, Inc.
20451 Seneca Meadows Parkway
Germantown, MD 20876-7005 USA
[eversensecgm.com](https://www.eversensecgm.com)

The Eversense® 365 Continuous Glucose Monitoring (CGM) System is indicated for continually measuring glucose levels for up to one year in people (18 years and older) with diabetes. The system is indicated for use to replace fingerstick blood glucose (BG) measurements for diabetes treatment decisions. Fingerstick BG measurements are required for calibration one time a week after day 13, and when symptoms do not match CGM information or when taking medications of the tetracycline class. The sensor insertion and removal procedures are performed by a health care provider. The Eversense 365 CGM System is a prescription device; patients should talk to their health care provider to learn more.

For safety information, see bit.ly/eversensesafety

Eversense, Eversense 365 Continuous Glucose Monitoring, and the Eversense logo are trademarks of Senseonics, Incorporated. Other brands and their products are trademarks or registered trademarks of their respective holders and are used solely for informative purposes. No relationship or endorsement should be inferred or implied. MKT-001723 Rev 2