

Share your CGM DATA

1 Allowing your healthcare provider access to your data

Clinic Share is an optional feature in the Eversense 365 App that enables you to securely share[†] your CGM data and reports with your healthcare provider's clinic by either sending an invitation or accepting a clinic request.

To invite your healthcare provider:


You will need their Eversense DMS ProClinic ID number. Please contact your HCP's office directly for this information.

Clinic ID: _____

1. Tap Menu > Share My Data > Clinic Share
2. Tap Clinic Requests > Send Requests
3. Enter your healthcare provider's clinic ID
4. Tap Add
5. Tap OK to confirm invitation

Once they accept the invitation, your healthcare provider will be able to view your data from their Eversense DMS Pro Account.



[†]The app applies appropriate safeguards to ensure your personal data is processed securely and in compliance with applicable laws.

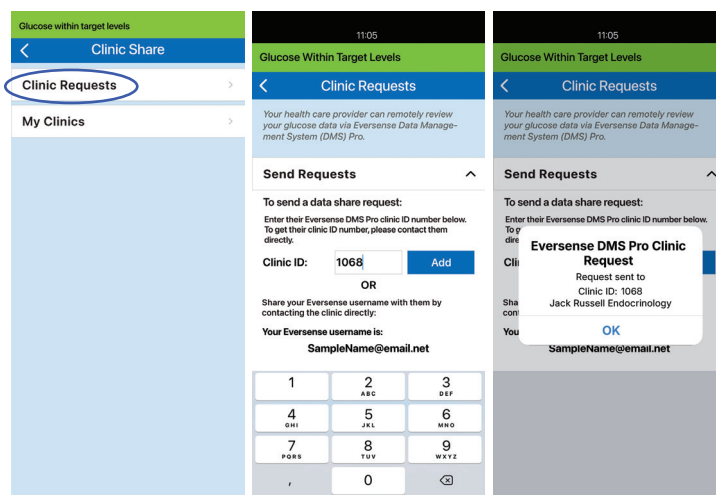

2 Sharing your data with family and friends

1. Tap Menu > Share My Data > My Circle
2. Tap Invite to My Circle > Add Email
3. Enter the email address for the person with whom you wish to share your data
4. Tap Send

They will need to download the Eversense NOW App to accept the invitation in the **Circles** menu to view your CGM data

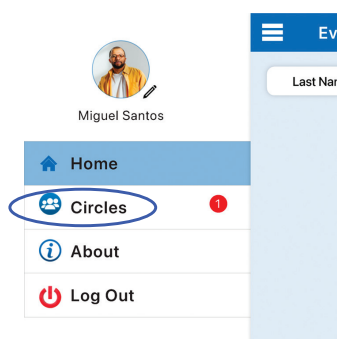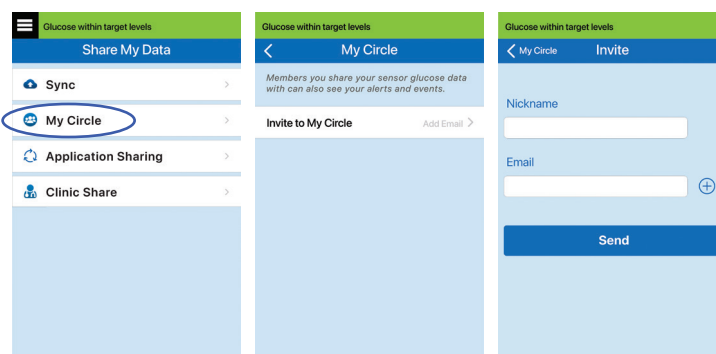

Scan this QR Code to watch a video on how to Share your CGM data with Healthcare Providers and Family



4 Accessing your CGM Data using a web browser


You can see your glucose data in a variety of reports and charts by logging into your Eversense account from a web browser on your computer, tablet, or mobile device.

1. Go to <https://us.eversensedms.com> and enter your log-in information.
2. On the left navigation pane, click on **Reports** to open the list of reports.

No need to remember multiple passwords! The mobile app and web browser use the same log-in information.

5 Viewing, Printing and Emailing Reports

Select **Reports** from the navigation bar on the left to view your reports. You may print or email the reports by clicking on the printer icon or the envelope icon in the right hand corner of the report.

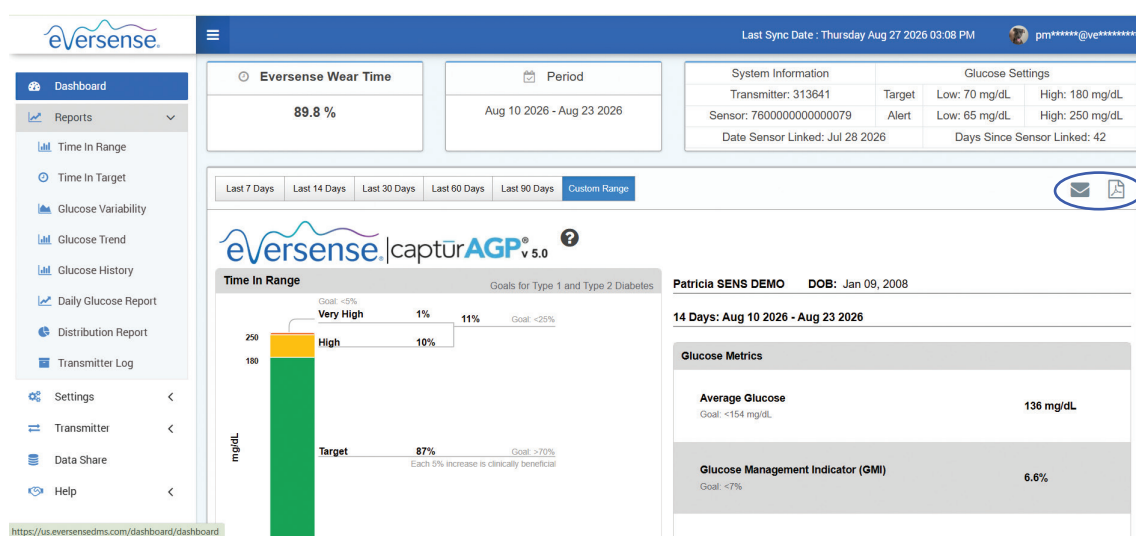

Located in the settings menu, REPORT SUMMARY settings allow you to set up a regularly-scheduled automated report that will be emailed to you. You can select how often the reports are sent, and which reports you want to receive.

To contact Customer Support, call 844-SENSE4U (736-7348)



© 2026 Senseonics, Inc.
All rights reserved.

Manufactured by Senseonics, Inc.
20451 Seneca Meadows Parkway
Germantown, MD 20876-7005 USA
eversensecg.com

The Eversense® 365 Continuous Glucose Monitoring (CGM) System is indicated for continually measuring glucose levels for up to one year in people (18 years and older) with diabetes. The system is indicated for use to replace fingerstick blood glucose (BG) measurements for diabetes treatment decisions. Fingerstick BG measurements are required for calibration one time a week after day 13, and when symptoms do not match CGM information or when taking medications of the tetracycline class. The sensor insertion and removal procedures are performed by a health care provider. The Eversense 365 CGM System is a prescription device; patients should talk to their health care provider to learn more.

For safety information, see bit.ly/eversensesafety

Eversense, Eversense 365 Continuous Glucose Monitoring, and the Eversense logo are trademarks of Senseonics, Incorporated. Other brands and their products are trademarks or registered trademarks of their respective holders and are used solely for informative purposes. No relationship or endorsement should be inferred or implied. MKT-002225 Rev 1